BEAUTIFLY

B-Lumi

Blush

series

Depilator laserowy

Laser hair remover

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User manual

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
2 OVERLAYS FOR DIFFERENT PARTS OF THE BODY INCLUDED
500,000 LIGHT PULSES
5 LEVELS
LARGE 5.2 CM² EPILATION ATTACHMENT WINDOW
EASY TO READ LCD DISPLAY
HANDY AND EASY TO USE
COSMETIC BAG INCLUDED

The B-Lumi cosmetic laser uses Intense Pulsed Light (IPL) technology, which is one of the most effective ways to quickly and easily remove unwanted hair and revitalize the skin. Compared to the usual method of laser hair removal, it is safer and more reliable, and the effects are more durable. Pulsed light hair removal requires the appropriate pulse wavelength. The B-Lumi has a light wavelength of 510nm-1100nm. The melanin in the hair absorbs high energy heat and the hair follicles naturally shrink and fall off. Using the B-Lumi device with Photon technology requires about 8 weeks for hair removal. After the first two treatments, you can notice hair reduction. After each epilation treatment, the hair thickness is reduced. The device can be used on most parts of the body.

Technical data	
Input power:	36W
Rated voltage:	100-240V
Rated frequency range:	50Hz - 60Hz
Storage temperature:	from 0°C to 40°C
Working temperature:	from 5°C to 35°C
Humidity at work:	from 30% to 85%

General specification

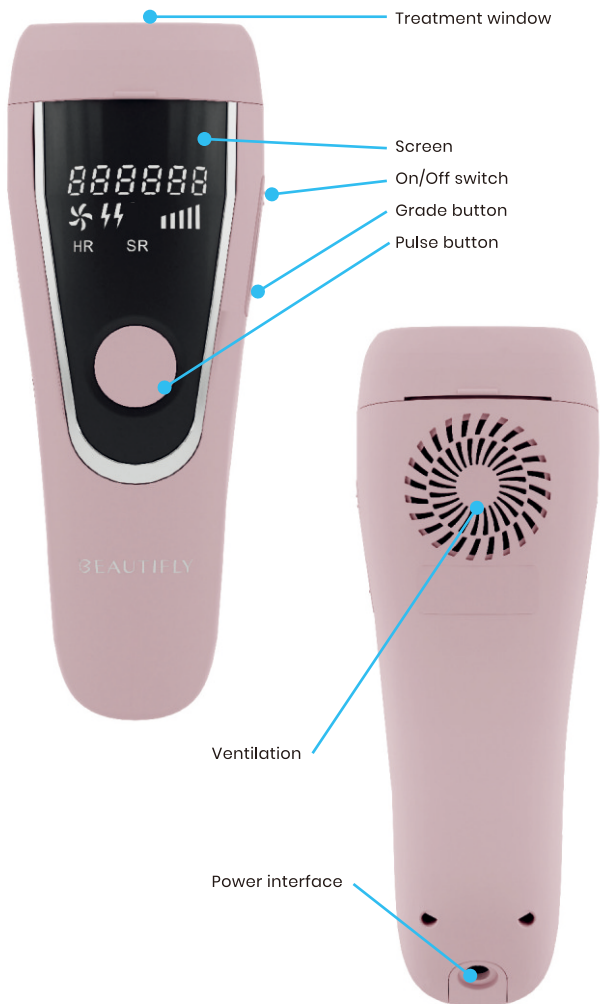
For the best results, we suggest that before the first use, read the instruction manual including all safety rules and warnings, for correct and safe using the device.

Function

B-Lumi has two functions and separate lamps for each of them.

- Function 1: IPL hair removal
- Function 2: IPL skin revitalization

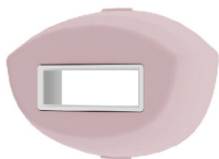Epilation

windows size: ~ 3.2 x 1.6 cm

Skin revitalization

windows size: ~ 3.2 x 1.9 cm

Automatic pulse mode

Press the pulse button for 2 seconds, the screen will show **⚡**, at the same time, the device will switch to automatic mode. The device will flash automatically when the pulse window is on your skin. When the machine is in auto mode, press the pulse button for 2 seconds, the device will switch to manual mode shown **⚡**.

Contraindications

1. The device is not intended for children under 15 years of age. Young people aged 15-18 can use the device with permission and / or under the supervision of parents or guardians.
2. Do not use if the skin in the area to be treated is tanned or burned. Using the device on tanned or burned skin can cause temporary or even permanent discoloration of the skin.
3. Do not use if you are pregnant, try to get pregnant, you suspect you may be pregnant or you are breastfeeding. The B-Lumi device has not been tested for this occasion, therefore, do not use it in these cases. Hormonal changes occurring during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
4. Do not use the device if you have a disease such as epilepsy or light-induced seizures or migraine. Energy flashes can cause severe headaches and migraines.
5. Do not use on areas covered with tattoos, dark brown or black places (such as dark freckles, birthmarks, or warts). Tattoos or dark spots can absorb too a lot of light energy that can cause swelling of the skin, temporary or permanent blisters, discoloration or whitening of the skin.
6. Do not use the device if you have a tendency to form keloids.
7. Do not use on the scalp, ears, neck, nipples or on mucous membranes (lips vulva smaller, around the vagina and anus). At those points of the bikini zone where the skin has stronger pigmentation, use a lower level of light intensity. Use in these places may cause swelling of the skin, blisters or temporary discoloration or skin bleaching.


8. Do not use if any active inflammation of the skin occurs in the treatment area, including acne, psoriasis, pimple, burns, herpes, abrasions or epidermis rash. These factors can make the skin more sensitive to light and make the procedure work unpleasant or cause swelling, blisters or infections on the skin. If the skin is enough healed, before using the device, consult a doctor.
9. Do not use the device if you are taking medicines that sensitize the skin to light including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazides, diuretics, sulfonyleureas, sulfonamides, DTICs, fluorouracil, vinblastine, griseofulvin, alpha-hydroxy acids (AHAS), Beta-Hydroxi Acids (BHAS), Retin-A, Accutane and / or local retinoids.
10. Do not use if you have a disease related to photosensitivity, e.g. porphyria, polymorphic light eruption, solar urticaria, lupus etc. These conditions make it the skin is more sensitive to light and the treatment can be unpleasant or causing swelling of the skin, blisters or infections.

Product principle

The B-Lumi hair removal device has integrated filters that protect the skin against harmful radiation and burn during the depilation procedure. Hair growth rate varies, depending on age, metabolism and other factors. Hair growth is usually divided into three stages.

1. **The growing phase** actively grows along with onions. At this stage, the concentration of melanin is the highest. Melanin determines the formation of hair color. Melanin concentration is also a decisive factor in the effectiveness of intensive pulsation light technology. Only hair in the growing phase will react to the liquid depilation process.
2. **Phase degeneration** with hair growth and rooting onion.
3. **A phase built** on the hair and fall off onions and fall out. Onion is in sleep until the next hair growth cycle.

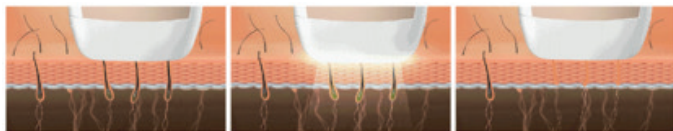The first stage

The second stage

The third stage

In the first stage of depilation, the B-Lumi device uses energy to act directly on the hair bulbs below the skin. Melanin in the middle and my hair bulbs absorb the impulse energy of light emitted by the device. The darker the color of the hair, the more light energy is absorbed, and this process stimulates the onion to enter the dormant phase.

In the second stage, during the depilation period, the hair will naturally fall out, and regeneration will be inhibited, which may take 1-2 weeks. Only hair in the growth phase will be removed. To effectively remove unwanted hair, repeat the depilation process once a week for a period of two months.

Operating warnings

1. The product will not turn on without inserting the lamp, insert the lamp before connecting to power. Do not block the intake air openings and the and ventilation holes during operation so that it does not adversely affect heat generation.
2. Wear protective goggles when using the device.
3. Clean the area of skin to be treated before use. It is not allowed to apply any cosmetics to avoid the destroy of instrument and affect result.
4. Adjust the energy level from low to high level. On condition of safety, the higher the energy level is, the better result it will be. After treatment, it is normal to have skin redness.
5. After the phototherapy treatment, use moisturizing preparation and sun block product.
6. After treatment then power off, please use wet cotton swab or cloth to clean treatment window. This can avoid hair and sebum of carbon sticking glass and affecting energy of shooting light.
7. Please don't let system and accessories touch liquid or put it in humid environment. Please store instrument in ventilated and dry environment and keep far away from children.
8. Please don't put metal object into heat emission hole. Otherwise there will be risk of electric shock.
9. Temperature control function: The device will automatically turn off the screen when the temperature is too high after using the device for a long time. You can turn on the screen again when it cooling down for 3 minutes.

Installation: use thumb and forefinger to hold the strip-like bulge of up and down sides, push into the slot to complete the installation.

Remove: use thumb and forefinger to hold the strip-like bulge of up and down sides, pull out along the slot.

Perform skin sensitivity test


If you are using the device for the first time or in a new area of the body, select a small area near the area you intend to undergo to test the sensitivity of the skin. Use the lowest level of the device. Follow the instructions for the procedure below, but emit no more than 2-3 flashes. After the skin test, wait 60 minutes and check that the skin does not react. If after 60 minutes the skin appears normal and there is no redness or other skin reaction, you are ready to perform the full procedure.

Choose the right intensity

Choose the suitable energy according to your skin tolerance and increase energy from one to five level gradually. If there isn't unsuitable or reaction, you could increase energy gradually.

Don't hold the applicator window still over one spot on the skin and flash multiple pulses on the skin. Doing so may cause too much heat and may result in skin swelling, blisters, or temporary or permanent changes in the color of your skin. After each flash, move treatment window slowly to a neighboring skin spot. Before returning to a previously treated skin spot, wait for at least 20 seconds.

Operation procedure

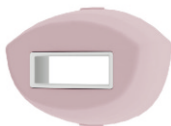Hair removal

(lamp depilation / detailed depilation)

Application

- Lip hair, armpit hair, body hair and legs, hair on places that affecting looks such as hair line on forehead and bikini area.
- The effect is different according to skin color. It is not suitable for white and light color hair.

Skin and hair color comparison table

skin color \ hair color	white	Ivory tint	beige	light brown	brown	dark brown and black
black	✓	✓	✓	✓	✓	✗
dark brown	✓	✓	✓	✓	✓	✗
brown	✓	✓	✓	✓	✓	✗
golden	✓	✓	✓	✓	✓	✗
red	✗	✗	✗	✗	✗	✗
white	✗	✗	✗	✗	✗	✗
gray	✗	✗	✗	✗	✗	✗

✓ applicable ✗ not applicable

How often to perform depilation treatment:

Usually, it is about 7–9 treatments. The first 3 treatments take place once a week. In the next 4–9 treatments, do one treatment every 2–3 weeks.

During the period of maintaining the depilation effect and preventing the hair roots from reactivating, regular treatments should be repeated so that the skin remains smooth. If the treatments are not continued, the hair roots will reactivate and the hair will grow back. In order to preserve the results, it is enough to repeat the treatments every 4–8 weeks, which will ensure that the hair roots will not activate and the hair regrowth will remain inhibited. The intervals between treatments may vary depending on the individual rate of hair regrowth.

How to perform depilation treatment:

- Step 1:** Before using the B-Lumi device, prepare the skin by shaving the selected area to be treated. The skin should also be washed.
- Step 2:** Put on safety glasses to avoid eye stimulation with IPL light.
- Step 3:** Select the desired intensity (see "Choose the right intensity") and start from a low level.
- Step 4:** Let the treatment window touches treated areas vertically, and treat the areas by turns. Do not repeat again and again in the same place, nor omit any area.
- Step 5:** Start in the initial area and repeat step 4, 2-3 times.
- Step 6:** After the treatment it is recommended to use a moisturizer after depilation.
- Step 7:** Turn off the device and use a damp cloth or swab to clean the treatment window.
- Step 8:** Place the device and accessories in a vanity case or box and store them in a safe place.

Photorejuvenation of the skin

(lamp wrinkle reduction)

Do you want to keep the skin young and healthy for years? We use our skin throughout the year by subjecting it to various atmospheric conditions. The exposure to the sun is particularly effective. Intense tanning usually leaves marks on the face or neckline, which is often difficult to mask with makeup. Similar discolorations also appear on the backs of the hands over time. How to deal with them?

People who want to visibly rejuvenate the skin and refresh its appearance, even out the color and reduce the stumbling discoloration, we offer photorejuvenation treatments with a modern B-Lumi cosmetic laser. Laser therapy B-Lumi allows for significant skin smoothing, complete removal of age spots or blishes.

The purpose of rejuvenation of selected parties are short and safe light pulses directed towards the surface of the skin. The light wave stimulates skin cells to produce collagen and elastin, which accelerate the regeneration of the skin. The effect of the therapeutic effect of IPL radiation is the improvement of blood supply and oxygenation of the skin, increase in the number of collagen and elastic fibers, better absorption of moisture, regulation of the dyes and keratinization process.

Effects:

- Elimination of signs of skin aging
- Reconstruction of collagen fibers
- Firm and elastic skin
- Reducing the pores
- Smoothing wrinkles
- Improving skin hydration
- Lightening discoloration

How often to perform photorejuvenation treatment:

To obtain the desired effects, a series of treatments should be carried out. It is recommended to carry out a series of 8 treatments, one treatment each week.

How to perform photorejuvenation treatment:

- Step 1:** Before proceeding with the treatment, thoroughly clean the skin with makeup and cosmetics. Also, wipe the lamp.
- Step 2:** Put on safety glasses to avoid eye stimulation with IPL light.
- Step 3:** Select the desired intensity (see "Choose the right intensity") and start from a low level.
- Step 4:** Allow the device window to touch selected areas vertically and treat them areas in turn. Do not repeat again and again in the same place, nor omit any area.
- Step 5:** Start in the initial area and repeat step 4, 2-3 times.
- Step 6:** After the procedure, it is recommended to use a strong moisturizer.
- Step 7:** Turn off the device and use a damp cloth or swab to clean the treatment window.
- Step 8:** Place the device and accessories in a case or box and store them in a safe place.