# BEAUTIFLY

Looky  
Young  
series

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
 /beautifuly.eu  @beautifuly.eu  beautifuly.eu

 /beautifulyeu  @beautifuly.eu

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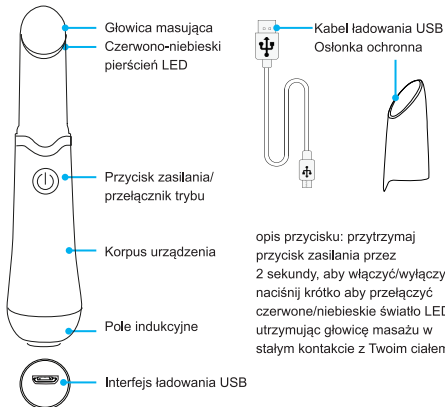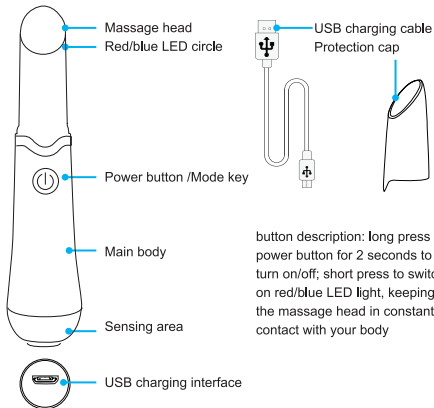
### Technical data

<b>Size:</b>	25*25*120mm
<b>Net weight:</b>	30g
<b>Temperature:</b>	42 ± 2°C
<b>Power:</b>	≤3W
<b>Charging voltage:</b>	== 5V
<b>Voltage:</b>	== 3.7V
<b>Charging current:</b>	160mA
<b>Battery capacity:</b>	200mA

## Product features

It integrates hot compress + high frequency micro-vibration + phototherapy + ion introduction. Perfect to be used with skin care products for efficient absorption of active ingredients to the bottom of the skin, ensuring best eye and facial beauty effects.


High-frequency vibration massage increases collagen, repairs elastic fibers, and reduces dark circles and bags under the eyes.

Magnetic thermal massage head with constant temperature hot compress ( $42^{\circ}\text{C} \pm 2^{\circ}\text{C}$ ), wakes up the skin and makes it firm, smooth and elastic. Blue light and red light therapy can penetrate into the bottom of the muscle effectively. It can enhance skin elasticity and dilute the melanin around the eyes.

Do not use the instrument while charging. The product will remind you when the battery is low. The red light blinks when the battery is low, and it continues to flash during charging. The solid blue light lights up when fully charged. When unused for 2 minutes, it goes into shutdown mode.

The modern design of the massage head adapts to the skin.

## Comparison diagram

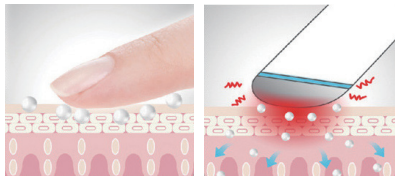
## How to use

Take an appropriate amount of skin care product and apply it around the eyes.

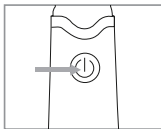According to the required skin care needs, long press for 2s to turn on the instrument and short press to switch between red and blue.

While massaging, hold the instrument so that the massage head is in contact with the skin and vibrating.

After using, dry the massage head with a tissue and cover it with a protective cover.

## How to use

**Eye care:** Apply eye cream and using the device slowly massage the skin from the inner to the outer corner of the eye. It will improve blood circulation and remove dark circles. Eye cream, that is deeply absorbed, can make the skin firm and supple.


**Forehead care:** Apply a skin care product and massage upward from the eyebrow to the hairline. Ion high-frequency vibration, deeply activates skin cell tissue. It stretches the skin above the nose, between the eyebrows, and on the forehead.

**Cheek care:** Apply a skin care product and massage the skin from the inside of the face outwards to promote blood circulation, achieve the effect of skin rejuvenation and dilute fine lines.

**Lip care:** Apply a lip balm along the contour of the lips. With a gentle massage, lift the corners of the mouth upwards to smooth the fine lines around them and improve hydration, making the lips tender and shiny.

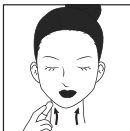**Neck care:** Repeat the massage with moisturizing cream from bottom to top. The ion effect and high frequency vibration massage are two major skin care magic weapons to revitalize the cells, deeply nourish the skin and dilute the neck lines.

## Storage

Keep this product at a dry place protected from humidity.

Keep this product out of reach of children.

## Safety instructions

Please read these instructions and follow them when using the device.

Do not use this product on wounded, bleeding, or swollen skin.

Do not massage eyeballs or inner parts of eyes.

Prevent this product from falling or strong collisions.

Do not place it in the shower room or humid areas.

Keep it away from strong magnet fields.

Do not expose it to sunlight.

Keep this product away from erosive fluids.

Never try to repair the product yourself. Repairs may only be performed by qualified specialists

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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The importer hereby declares that the Beautifly Looky Young device complies with the Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, and with the Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the Declaration of Conformity is available at: (link below)

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