
medivon.
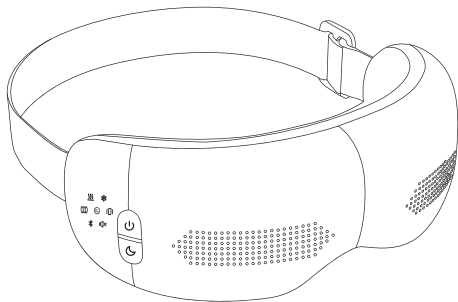
# User manual

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## Device Manual

Fully charge the device before using it for the first time and after a long period of non-use. After charging the battery, you can hold the ON/OFF button for at least 2 seconds to start the massager.

When the massager is turned on, pressing the power button cyclically changes the day massage modes, and the moon button cyclically switches between night modes.

The diodes on the device indicate which functions are currently active in a given mode: cooling, heating, vibration.

A quick double press of the moon button plays the next track, while holding it down mutes the music.

Recharge the device when the battery level is low. Do not fully discharge the battery. The supply voltage is 5V/1A. Do not use the device while it is charging.

After finishing the massage, the device can be turned off by pressing the On/OFF button again.

To connect to the massager and play your own music, just search for the device in the list of visible bluetooth devices and connect. The music will play automatically when you select a song.

The strength of the massage should be adjusted to your own preferences. The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use. Pregnant women cannot use the massager.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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The manufacturer hereby declares that the Horizon Lynx device complies with Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address:  
[www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska  
Model: Horizon Lynx