

# Panasonic


Operating Instructions  
(Household) IPL Hair Removal System  
Model No. ES-WH90

CE

English.....	5	Español .....	85
Deutsch .....	25	Polski.....	105
Français.....	45	Türkçe.....	125
Italiano .....	65		
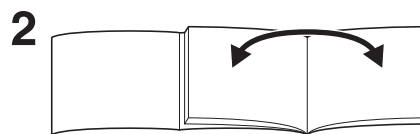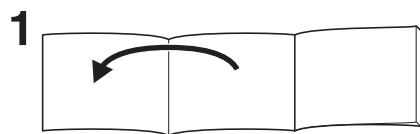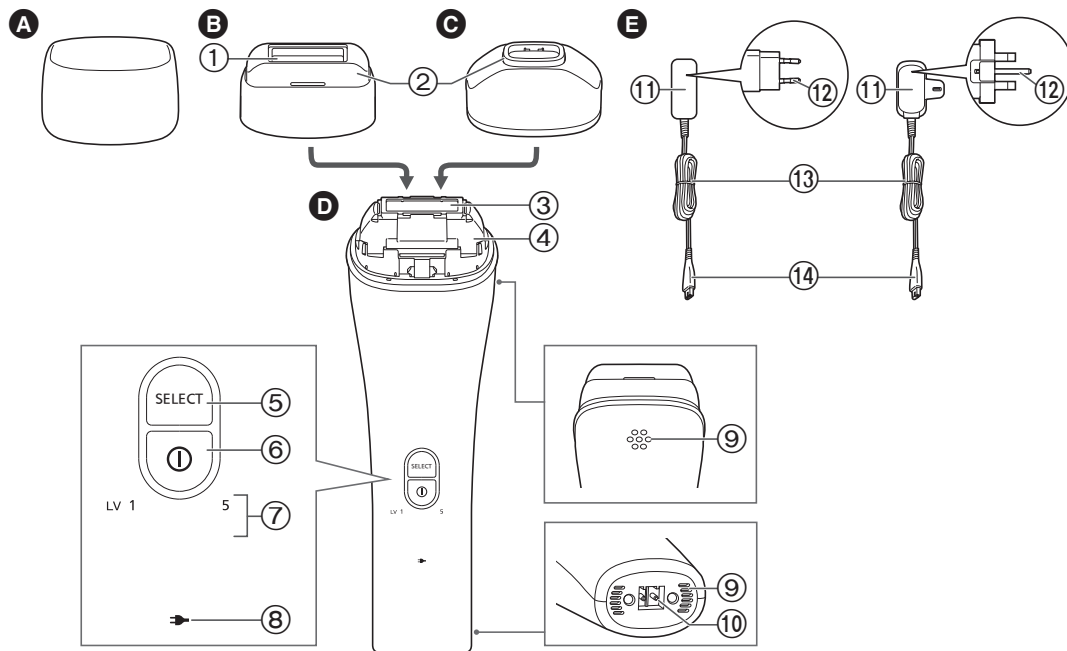
# Panasonic

Operating Instructions  
(Household) IPL Hair Removal System  
Model No. ES-WH90

## Contents

<b>Safety precautions.....</b>	<b>8</b>	<b>For healthy skin .....</b>	<b>18</b>
<b>Intended use .....</b>	<b>12</b>	<b>Lamp indication .....</b>	<b>18</b>
<b>Parts identification .....</b>	<b>13</b>	<b>Cleaning .....</b>	<b>19</b>
<b>Treatment areas .....</b>	<b>13</b>	<b>Q&amp;A.....</b>	<b>20</b>
<b>Guide for treatment frequency, time ...</b>	<b>14</b>	<b>Troubleshooting .....</b>	<b>21</b>
<b>Charging.....</b>	<b>14</b>	<b>Battery life .....</b>	<b>22</b>
<b>Before use .....</b>	<b>15</b>	<b>Removing the built-in rechargeable batteries.....</b>	<b>22</b>
<b>When using the appliance for the first time .....</b>	<b>16</b>	<b>Specifications .....</b>	<b>23</b>
<b>How to use .....</b>	<b>16</b>		

Thank you for purchasing this Panasonic product.

**Before operating this unit, please read these instructions completely and save them for future use.**

## Warning

- This appliance can be used by children aged from 18 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The supply cord cannot be replaced. If the cord is damaged, the AC adaptor should be scrapped.
- Keep the appliance dry.
- Do not use anything other than the supplied AC adaptor.

- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

## Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

### Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

 <b>DANGER</b>	Denotes a potential hazard that will result in serious injury or death.
 <b>WARNING</b>	Denotes a potential hazard that could result in serious injury or death.
 <b>CAUTION</b>	Denotes a hazard that could result in minor injury or property damage.

The following symbols are used to classify and describe the type of instructions to be observed.

 This symbol is used to alert users to a specific operating procedure that must not be performed.
 This symbol is used to alert users to a specific operating procedure that must be followed in order to operate the unit safely.

## **WARNING**

### ► To prevent skin or health problems

**Persons using an implanted medical electronic device, such as a pacemaker, should not allow this unit to touch the skin where the medical electronic device is implanted.**

- Doing so may cause an accident or illness.

**Do not use during pregnancy or menstruation.**

- Doing so may cause skin problems due to unstable hormone balance.

**Do not use on tattoos.**

- Doing so may cause keloids by burns or redness or other rashes on skin.

**Do not shine directly into the eyes.**

- Doing so may cause injury to the eyes.

**Only use on underarms, arms, legs, bikini-line or face (when using the face attachment).**

- Failure to do so may cause burns or redness or other rashes on skin.

### ► Power supply

**Do not connect or disconnect the adaptor to a household outlet with a wet hand.**

- Doing so may cause electric shock or injury.

# WARNING

**Do not damage or modify, or forcefully bend, pull, or twist the cord. Also, do not place anything heavy on or pinch the cord.**

 **Never use the appliance if the AC adaptor is damaged or if the power plug fits loosely in a household outlet.**  
- Doing so may cause electric shock or fire due to a short circuit.

 **Do not use anything other than the supplied AC adaptor. Also, do not charge any other product with the supplied AC adaptor.**  
- Doing so may cause burn or fire due to a short circuit.

 **Do not use in a way exceeding the rating of the household outlet or the wiring.**  
- Exceeding the rating by connecting too many plugs to one household outlet may cause fire due to overheating.

 **Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the AC adaptor.**

**Fully insert the adaptor or appliance plug.**  
- Failure to do so may cause fire or electric shock.

 **Always unplug the adaptor from a household outlet when cleaning.**  
- Failure to do so may cause electric shock or injury.

 **Regularly clean the power plug and the appliance plug to prevent dust from accumulating.**  
- Failure to do so may cause fire due to insulation failure caused by humidity.  
Disconnect the adaptor and wipe with a dry cloth.

## ► This product

 **This product has a built-in rechargeable battery. Do not throw into flame or heat. Do not charge, use, or leave in high temperatures.**  
- Doing so may cause overheating, ignition, or explosion.

 **Do not store within the reach of children. Do not let them use it.**  
- Doing so may cause an accident or injury.

 **Do not modify or repair.**  
- Doing so may cause fire, electric shock, or injury.  
Contact an authorized service centre for repair (battery change etc).

 **Never disassemble except when disposing of the product.**  
- Doing so may cause fire, electric shock, or injury.

 **Do not use if the air inlet is blocked.**  
- Doing so may cause fire or burns.

## ► Main body and AC adaptor

 **Do not store the main body and AC adaptor in a bathroom or location with high humidity. Also, do not place in a location where water may spill (such as top of a vanity).**  
- Doing so may cause electric shock or fire.

 **Do not immerse the main body and AC adaptor in water or wash it with water.**  
 **Do not use the main body and AC adaptor in the bathroom or location with high humidity.**  
- Doing so may cause electric shock or fire due to a short circuit.

GB

D

F

I

E

PL

TR

## **WARNING**

### ► In case of an abnormality or malfunction

**Immediately stop using and remove the adaptor if there is an abnormality or malfunction.**

- Using it in such conditions may cause fire, electric shock, or injury.

**<Abnormality or malfunction cases>**

- **The main unit, adaptor or cord is deformed or abnormally hot.**
  - **The main unit, adaptor or cord smells of burning.**
  - **There is abnormal sound during use or charging of the main unit, adaptor or cord.**
- Immediately request inspection or repair at an authorized service centre.

## **CAUTION**

### ► To protect the skin

**The following persons should not use this appliance:**

- Persons with dermatitis
- Persons suffering from atopic dermatitis or rashes and other skin diseases
- Persons with allergic diathesis, or sensitive skin which is easily irritated by cosmetics, clothing, metals, etc.
- Persons with skin that is sensitive to sunlight

**Do not use in the following body area:**

- On suntanned skin
- Dark brown or black skin areas such as moles (Refer to the sheet for "Suitable Hair and Skin Tones".)
- Scabs, injuries, rashes, bruises, warts, pimples, skin cancer, ringworm (dermatophytosis), folliculitis, irritated skin, insect bites, etc.
- Angioma
- Areas with cosmetic surgery

**Do not use in the following cases:**

- When you are not well or when you feel an abnormality in your skin or body

**Do not use on the same area repeatedly. (Try using the appliance up to once every 3 days even if removal of unwanted hair is not immediately obvious.)**

- Doing so may cause injury to the skin, inflammation, bleeding, increased pain or other problems or aggravated symptoms.

**Do not remove the attachment except during cleaning.**

- Doing so may cause injury to the skin or fingers.

# CAUTION

 **Only for use by persons over the age of 18.**  
- Failure to do so may cause injury or inflammation of the skin, bleeding or increased pain.

 **Before use, check that there is no distortion, cracks or damage to the light guide and flash window.**  
- Failure to do so may cause injury to the skin.

 **Persons on medication should consult a dermatologist.**  
- Skin may become sensitive to light stimulation depending on the type of medication being taken, and may cause skin problems.

 **If inflammation (such as skin turning red) occurs, stop using the appliance, and consult a dermatologist if the inflammation continues even after 2 days.**  
- Failure to do so may cause symptoms to worsen.

 **Use after removing metal objects such as accessories.**  
- Failure to do so may cause an accident or injury.

## ► Power supply

 **Do not allow metal objects or trash to adhere to the power plug or the appliance plug.**  
- Doing so may cause electric shock or fire due to a short circuit.

 **Do not wrap the cord around the adaptor when storing.**  
- Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.

 **Disconnect the adaptor or the appliance plug by holding onto the adaptor or the appliance plug instead of the cord.**  
- Disconnecting by holding onto the cord may cause electric shock or injury.

 **Disconnect the adaptor from the household outlet when not charging.**  
- Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

## ► This product

 **Do not use on wet skin.**  
- Doing so may cause an electric shock.

 **Do not drop or subject to shock.**  
- Doing so may cause injury.

 **Do not remove the attachment or touch the flash window immediately after using IPL hair removal.**  
- Doing so may cause burns as the flash window becomes hot.

## ► Note the following precautions

 **Do not use on clothes or head hair.**  
- Doing so may cause burns or discoloration.

GB

D

F

I

E

PL

TR

## Handling of the removed battery when disposing

### **DANGER**

The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products. Do not charge the battery after it has been removed from the product.

- Do not throw into fire or apply heat.
  - Do not hit, disassemble, modify, or puncture with a nail.
  - Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
  - Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
  - Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
  - Never peel off the tube.
- Doing so may cause overheating, ignition, or explosion.

### **WARNING**

After removing the rechargeable battery, do not keep it within reach of children and infants.

- The battery will harm the body if accidentally swallowed. If this happens, consult a physician immediately.

**If the battery fluid leaks out, do not touch the battery with your bare hands.**

- The battery fluid may cause blindness if it comes in contact with your eyes.

Do not rub your eyes. Wash immediately with clean water and consult a physician.

- The battery fluid may cause inflammation or injury if it comes in contact with the skin or clothes.

Wash it off thoroughly with clean water and consult a physician.

### **Intended use**

- Please observe the following before performing IPL hair removal.
  - This appliance should only be used on bare skin with no cosmetics, lotions etc.

When using the lotion or shaving agents (gels, etc.), rinse off the agent thoroughly with water and let the skin dry.

· Be sure to shave the area (body or face).

Before use, wipe away any water or sweat from the treatment area, and any hair sticking to the skin first.

Hairs or shaving agents stuck to the skin can foul the light guide. It may become difficult to clean the light guide.

- Do not use alcohol, nail polish remover, or detergent (ex. liquid hand soap), etc. to clean the appliance.

Doing so may cause failure, cracking or discoloration of the main body. Wipe with a dry cloth.

- Do not leave the appliance in an area with high temperature. Doing so may cause malfunction.


## Parts identification

- A** Cap
- B** Body attachment
- C** Face attachment
  - ① Light guide
  - ② Flash switch
- D** Main body
  - ③ Flash window
  - ④ Flash lamp
  - ⑤ Level selector switch (flashes)
  - ⑥ Power switch
  - ⑦ Level indication lamp (blue) (flashes)  
(See pages 16 and 18.)
- ⑧ Charging indication/ notification lamp (red)  
(See pages 14 and 18.)
- ⑨ Air inlet
- ⑩ Socket
- E** AC adaptor (RE7-87)  
(The shape of the AC adaptor differs depending on the area.)
  - ⑪ Adaptor
  - ⑫ Power plug
  - ⑬ Cord
  - ⑭ Appliance plug

## Removing and mounting the body/face attachment

### ► Removing the attachment

Holding the attachment by side, remove from the main body.
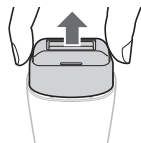
### ► Mounting the attachment

Push the attachment until it clicks.

## Treatment areas


### Face

Do not use on the  area, such as around the eye, eyebrow and lip, etc.

- Tie up hair or fix it with a hair band or similar. Also, please also fix short hair with a clip or similar.
- Please take care in particular of sunburn before and after using on your face.
- Do not use for the man's beard.

### Underarms

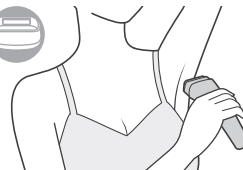
### Arms

### Bikini-line \*1,2

### Legs

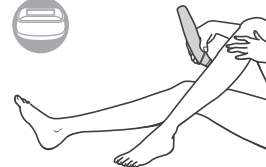\*1 For bikini-line treatment, wear a swim wear, under wear etc. first to check the treatment area.

\*2 Do not use in the abdomen and genitals.

GB

D

F

I

E

PL

TR

## Guide for treatment frequency, time

### Frequency

#### Legs/Arms

Use once every 2 weeks for the first 2 months.

Use once every 4 weeks from the 3rd month onward.

#### Underarms/Bikini-line/Face

Use once every week for the first 2 months.

Use once every 2 weeks from the 3rd month onward.

### Guide for treatment time

Legs	Approximately 10 minutes per session for both legs (approximately 190 flashes)
Arms	Approximately 5 minutes per session for both arms (approximately 90 flashes)
Underarms	Approximately 2 minutes per session for both underarms (approximately 30 flashes)
Bikini-line	Approximately 2 minutes per session (approximately 30 flashes)
Face	Approximately 3 minutes per session (approximately 40 flashes)


- Note that, for the facial brightening, it is recommended to use the appliance once every 3 days.
- Do not use on the same area repeatedly, even if removal of unwanted hair is not immediately obvious. Try using the appliance up to once every 3 days. (Treatment effects may differ depending on the area or individual body types.)
- If the appliance has not been used for a while, resume use at the same frequency as the first 2 months.
- Treatment time may differ depending on the area or amount of hair, treatment area and how the appliance is used.

14

- Please shave any unwanted hair first from the treatment area. Please use a separate razor or electric shaver.
- The appliance differs from permanent hair removal treatment. Unwanted hair will gradually grow back if treatment is stopped.

## Charging

- 1 Insert the appliance plug firmly into the socket at the bottom of the main body.**
- 2 Firmly insert the adaptor into a household outlet.**

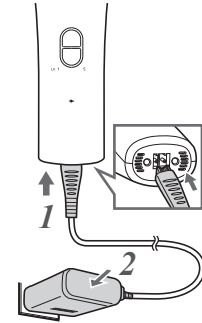- 3 Check that the charging indication/notification lamp glows.**

While charging	After charging is completed
 Red light glows.	 Red light goes off.

- Charging is completed in about 3 hours.
- The charging time may differ with the charging capacity.

<To check if charging is completed>

Press the power switch with the AC adaptor connected. The charging indication/notification lamp glows, and goes off after approximately 5 seconds.

#### **4 Remove the adaptor from the household outlet, and remove the appliance plug from the main body.**

(for safety and to reduce energy usage)

- There will be no problems with the main body even if it is charged for 3 hours or more.

#### **Notes**

- The charging time may differ when the appliance is being used for the first time, or if it has not been used for half a year or more. The lamp on the main body may not glow for several minutes after charging starts, but will glow if charging continues.
- The appliance can be used for approximately 600 flashes with a full single charge.
- The number of flashes may differ depending on the ambient temperature and how the appliance is used.
- Recommended ambient temperature for charging is 15 – 35 °C. Charging may not be possible outside of the recommended temperature.
- If there is crackling noise on radios or other devices while charging, charge the appliance using a different household outlet.
- Charging the appliance every time it is used will not affect the life of the rechargeable battery.

#### **Before use**

- Do not use the appliance in dark places. Please close your eyes when flashing your face with the face attachment. Doing so may cause temporary irritation in the eyes.
- This appliance should only be used on bare skin with no cosmetics, lotions etc. When using the lotion or shaving agents (gels, etc.), rinse off the agent thoroughly with water and let the skin dry.
- Be sure to shave the area (body or face). Before use, wipe away any water or sweat from the treatment area, and any hair sticking to the skin first. Hairs or shaving agents stuck to the skin can foul the light guide. It may become difficult to clean the light guide.

## When using the appliance for the first time

Check that the output level suits your skin.  
Treated skin may temporarily turn red or feel painful depending on the body type or condition of the person using the appliance.

**Set the output level to 1 (LV1) and flash the treatment area once.**

After 2 days ↓ If there are no skin abnormalities


**Set the output level to 1 (LV1) and start treatment on the area that you want.**

• Guide for treatment frequency, time (See page 14.)

After 2 weeks\* ↓ If there are no skin abnormalities  
\*For armpits/bikini line/face: after 1 week

**To increase the output level for treatment, flash once on the treatment area, and check your skin condition for around 2 days before further treatment. Increase the output level one step at a time.**


• If there is intense pain or redness after 2 or more days, stop using the appliance and consult a dermatologist.
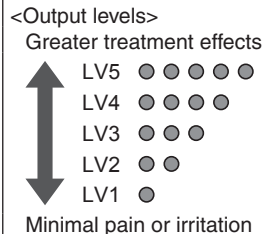
## How to use

- The power will not turn on even when pressing the switch with the AC adaptor connected to the main body.
- The appliance may not operate in an ambient temperature approximately 5 °C or lower.


**1 Remove the cap.**

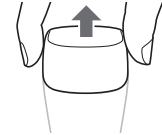**2 Attach the attachment to suit the area where you want to use it.**

**3 Press and hold (for approximately 3 seconds) to turn the power on.**

- The power will turn off automatically if the appliance is not used for approximately 1 minute.

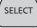**4 Press  to select the output level.**

- Level 1 is selected when the power is turned on.
- The level cycles through 1 → 2 → 3 → 4 → 5 → 1...each time the level selector switch is pressed.
- The level indication lamp blinks, and then turns on. (A beep sounds, indicating that preparation for flash is complete.)


## 5 When the level indication lamp glows, apply the flash switch closely against the skin and keep it pressed until flashes.

When the appliance beeps, press the head perpendicularly against the skin.

- Please do not press too strongly.
- Use the appliance after checking that there are no other people nearby.
- Light might leak from the air inlet on the back of the main body when flashing the treatment area; however, there will be no problem.
- Points for IPL hair removal (See this page.)
- When the appliance beeps, preparation for flash is complete.

## 6 Press to turn the power off.


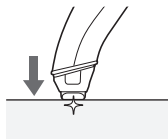
### Points for proper IPL hair removal

#### Point 1


**Apply the flash switch perpendicularly against the skin and keep it pressed against the skin.**

The flash window will not flash unless it is applied perpendicularly to the skin.
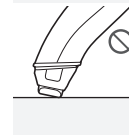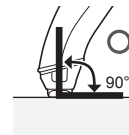
#### Point 2


The size of the  area is a guide for the IPL hair removal range of each flash.

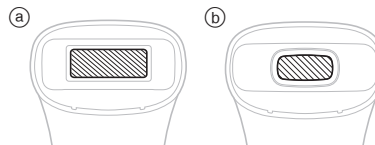(a) Body attachment

(b) Face attachment

- When the appliance beeps, preparation for flash is complete.

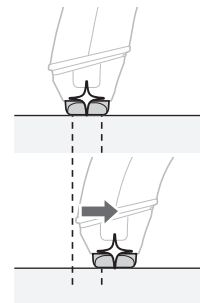GB

D

F

I


E

PL

TR

## Continuous flashes

After one flash, slide the appliance to the next flash area while it is still pressed against the skin. Up to 10 continuous flashes at an interval of approximately 3 seconds is possible.

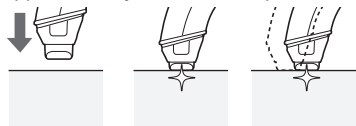- The flash window will no longer flash after 10 flashes. Remove the appliance from the skin and try again.
- Do not use the appliance continuously on the same area.
- Depending on usage conditions, the flash interval may vary.

## For healthy skin

The following care is recommended after IPL hair removal.

1. Cool the skin using a towel soaked in ice water.

2. Moisturize the skin with the lotion etc.

## Lamp indication

### ► Normal indication

	<p>The level indication lamp (blue) blinks (once/second). The image on the left is when level 5 is selected.</p> <p>➔ Preparation for flash in progress.</p>
	<p>The level indication lamp (blue) glows. The image on the left is when level 5 is selected.</p> <p>➔ Flashing is possible.</p>

### ► Charging notification

	<p>LV 1                      5</p> <p>The charging indication/notification lamp (red) blinks (once/second).</p> <p>The battery level is low. Charging is recommended.</p>
--	---

### ► Malfunction display

	<p>LV 1                      5</p> <p><b>[Abnormal charging]</b> The charging indication/notification lamp (red) blinks quickly (twice/second).</p> <p>Recommended ambient temperature for charging is 15 – 35 °C. Check the room temperature of the place where the appliance is going to be charged, and charge it again.</p>
--	---

### [Other malfunctions]

The level indication lamp (blue) and charging indication/notification lamp (red) blink simultaneously (once/second), the appliance beeps twice, and the power turns off after 5 seconds.

The appliance may not operate in an ambient temperature approximately 5 °C or lower. Check the room temperature of the place that the appliance is going to be used, and turn it on again.

(The appliance may not adapt to the room temperature immediately after changing the location.)

If the symptom occurs repeatedly, contact an authorized service centre.


## Cleaning

Always turn the main body power off and wait for 2 to 3 minutes before cleaning.

### ► Flash window/Light guide

1. Remove the attachment.

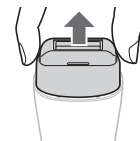2. Wipe any dirt clean carefully with a dry cloth.

Flash window

Light guide



3. Attach the attachment.

Insert until it clicks.

### ► Air inlet

Wipe any dirt clean with a dry cloth.

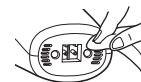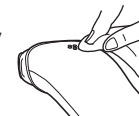GB

D

F

I

E

PL

TR

Q&A	
Question	Answer
What type of light is emitted from the flash window?	▶ Instantaneous flash from a xenon lamp.
Why can the appliance be used only by those over the age of 18?	▶ This is because of unstable hormone balance when growing up and the skin can become sensitive easily.
On which areas can the appliance be used?	▶ Only use the appliance on underarms, arms, legs, bikini-line and face (when using face attachment). (Do not use the body attachment on the face or other areas.)
Is it alright to use the appliance everyday?	▶ Using the appliance on the same area everyday is not recommended. Please follow the guide for treatment frequency and time. (See page 14.)
Why is there a burning smell after a flash?	▶ This is because the flash applies heat to unwanted hair, and is not a sign of a malfunction. Shave any hair if it has grown too long. Also remove all trimmings remaining on the skin.
The light feels very bright while using the appliance but is this a problem?	▶ Apply the flash switch closely against the skin. (See page 17.) If the light is still a concern, wear sunglasses or other eye protection.

20

Question	Answer
Can the flash lamp be replaced?	▶ The flash lamp cannot be replaced. If the flash lamp is damaged or the like, please contact the store where you purchased the unit or a service centre authorized by Panasonic for repair.
The light feels hot while using the appliance but is this a problem?	▶ The sensation is different for each person but it may feel hotter if hair has not been shaved properly. ▶ Treated skin may temporarily turn red or feel painful depending on the body type or condition of the person using the appliance. Stop using the appliance immediately, and keep the skin cool and moist.
The treated skin turns red while using the appliance but is this a problem?	▶ If the symptoms subside within 2 days, use the appliance at a level that the symptoms do not appear. If the symptoms persist after 2 days, consult a dermatologist. (See page 16.)
The pain and redness do not subside after use but is this a problem?	▶ Regardless of the remaining capacity of the battery, the flash power is the same.
When the remaining capacity of the battery is low (charging indication/notification lamp blinks), does the flash power drop?	

Question	Answer
Can the appliance be used by men?	▶ Men can use the appliance apart from on the beard, but because generally men have more and thicker body hair than women, the stimulation of the skin may be stronger.
Dirt on the light guide cannot be removed.	▶ If the dirt cannot be removed even after cleaning, contact the store where you purchased the unit or a service centre authorized by Panasonic for repair. [Repair part] Body attachment: WESWH90W3108

## Troubleshooting

Problem	Action
The flash window does not flash.	▶ Check that the level indication lamp (blue) has changed from blinking to glowing. ▶ Press the flash switch perpendicularly against the skin until it flashes. (See page 17.) ▶ If 10 continuous flashes have been emitted, temporarily move the appliance away from the skin.
The flash preparation time is long.	▶ The flash preparation time may become longer to prevent the main body from overheating.
The power turned off after all lamps blinked and the appliance beeped twice.	▶ The appliance may not operate in an ambient temperature approximately 5 °C or lower. (See page 19.)
The power turned off during use.	▶ For safety reasons, the power will turn off automatically if the appliance is not used for approximately 1 minute after the power is turned on. Turn the power on again.
The appliance can only be used a few times after charging.	▶ The battery has reached the end of its life. (See page 22.)

CB

D

F

I

E

PL

TR

Problem	Action
The appliance does not operate even after pressing the power switch.	▶ Press and hold the power switch (for approximately 3 seconds) when turning the power on.
	▶ Charge the appliance. (See page 14.)

If the problems still cannot be solved, contact the store where you purchased the unit or a service centre authorized by Panasonic for repair.

## Battery life

Approximately 3 years.

If the operating time is significantly shorter even after a full charge, the battery has reached the end of its life.

(The battery life may change due to usage and storage conditions.)

The battery in this appliance is not intended to be replaced by consumers. Have the battery replaced by an authorized service centre.


## Removing the built-in rechargeable batteries

### Remove the built-in rechargeable batteries before disposing of the appliance.

Please make sure that the batteries are disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it may cause malfunction.

- Remove the appliance from the AC adaptor.
- If the appliance can be turned on and emit light, use it until the battery is depleted.
- Wait for at least 10 minutes after turning the power off before disassembling the appliance.
- Perform steps ① to ⑦ and lift the battery, and then remove it.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.
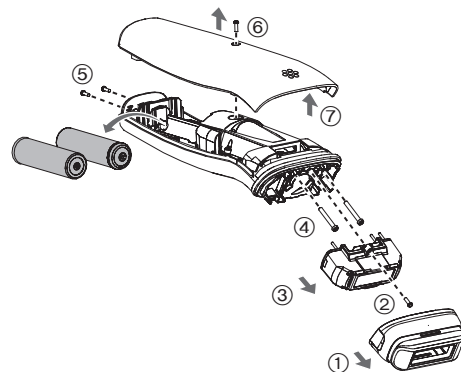
### For environmental protection and recycling of materials

This appliance contains Li-ion batteries.  
Please make sure that the batteries are disposed of at an officially assigned location, if there is one in your country.

### Specifications

Power supply	See the name plate on the AC adaptor. (Automatic voltage conversion)
Charging time	Approx. 3 hours
Airborne Acoustical Noise	Output level 5: 47 (dB (A) re 1 pW)

This product is intended for household use only.

### Disposal of Old Equipment and Batteries Only for European Union and countries with recycling systems

These symbols on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries must not be mixed with general household waste.

For proper treatment, recovery and recycling of old products and batteries, please take them to applicable collection points in accordance with your national legislation.

By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment.

For more information about collection and recycling, please contact your local municipality.

Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

#### Note for the battery symbol (bottom symbol)

This symbol might be used in combination with a chemical symbol. In this case it complies with the requirement set by the Directive for the chemical involved.

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© Panasonic Corporation 2016
Panasonic Corporation