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MASAŻER DO STÓP MMS-01

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
IMPORTANT SAFEGUARDS

- Please read all instructions carefully before using your appliance.
- Take special care when using near children.
- Do not let the power cord hang over sharp edges or touch hot surfaces.
- Do not use the appliance if the power cord or plug is damaged. In such case, return the appliance to an authorised service centre.
- Using attachments not recommended by the manufacturer may result in damage to the appliance, fire or personal injury.
- Do not immerse the power cord, plug or main body in water or any other liquid.
- Do not use appliance in close vicinity of bathtub or shower. Massager should be placed at least 1 m away from them.
- Always switch off device before unplugging it.
- Never leave the appliance unattended when in use.
- The appliance must only be used for its intended purpose.
- The water level should not exceed the MAX mark on the inner surface of the device.
- Never operate the device with wet hands.
- Do not stand up in the unit, use it only when you are sitting.
- Use only tap water and non-foaming bath salts and essences as the foam can clog the device.
- Do not use the device on soft surfaces, e.g. carpets as this may block the air flow.
- Do not move or unplug device during use.
- After use, remove water using the water spout at the back of the device.
- Special care must be taken during use of a device by the people non-sensitive for temperature.
- This appliance is not intended for use by children under 8 years of age and individuals with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved. Cleaning and user maintenance shall not be carried out by children without supervision.
- Store the appliance and its cord out of reach of children under 8 years of age.
- Make sure that children do not play with the appliance.
- For home use only. Do not use for medical purposes.
- Always unplug the appliance from the power supply after use and before cleaning.
- Do not plug the power cord with wet hands.
- This appliance is intended for household use only.
- Do not remove the plug from the mains by pulling by the cord.

- ▶ Do not use the appliance outdoors.
- ▶ Keep your children safe: do not leave any parts of the packaging readily accessible (i.e. plastic bags, cardboard boxes, polystyrene inserts, etc.).
- ▶ **WARNING! Never let children play with plastic film or bags. Choking/asphyxiation hazard!**

FEATURES

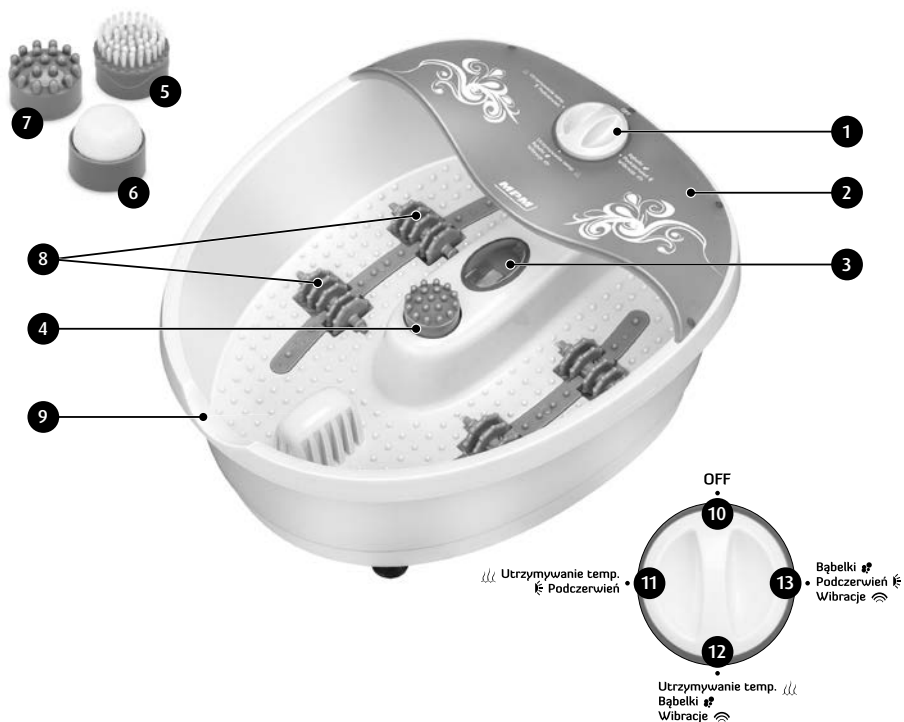1. ON/OFF and operation mode switch
2. Anti-splash lid
3. Infra-red radiator
4. Mounting hole for massage tips 5 6 7

5. Brush
6. Pumice
7. Massage ball
8. 4 removable foot massage rollers
9. Water spout

OPERATION MODES

10. OFF
11. Warm/Infrared
12. Warm/Bubbles/Vibration
13. Bubbles/Infrared/Vibration

USING THE APPLIANCE

- ▶ Appliance can be used for dry vibration massage or wet massage.
- ▶ Dry vibration massage can improve blood circulation and release muscle tension.
- ▶ Wet massage (with optional addition of non-foaming bath salts) helps to calm and relax.
- ▶ Infra-red stimulates blood circulation.

WARNING! Do not use the massager if you are pregnant, or suffer from thrombosis, open wounds, broken skin and varicose veins. Consult your doctor first!

1. Make sure that device is unplugged.
2. Fill the massager with warm water.
3. The water level should not exceed the MAX mark on the inner surface of the device.
4. Plug the device into the socket.
5. Place the massager on even surface. Place the feet inside the massage tub and use the switch to select the operating mode.

MESSAGE IN WARM WATER WITH INFRARED HEATING

– the built-in heating element will keep the water warm during the massage session.

VIBRATION MESSAGE WITH BUBBLES IN WARM WATER

– the vibration, bubbles and the air coming through the outlets will generate additional relaxing effect and improve blood circulation.

VIBRATION MESSAGE WITH BUBBLES AND INFRARED HEATING

– the vibration, bubbles and infra-red will generate additional relaxing effect and improve blood circulation, especially in long winter evenings.

MESSAGE TIPS

Massager is equipped with:

- ▶ 3 massage tips **5** **6** **7**, to be fitted in the mounting hole **4** in the middle part of the massager.
- ▶ 4 removable massage rollers **8**. Massage can be performed with or without the rollers.

MESSAGE DURATION

1. 10 minutes of a foot massage a day is enough to release psycho-physical stress.
2. When needed, you can change duration and frequency of a massage.

WARNING! The massager is designed for intermittent use. It can be applied without interruption not longer than 20 minutes, then the motor must be allowed to cool.

CLEAN AND CARE

1. Switch the appliance off and ensure that the massager has completely cooled. Unplug it from wall socket.

2. Remove water using water spout 9.
3. Rinse out with water. You can use water with a mild detergent, especially after using bath salts. Do not use hard abrasives or cleaners.
4. Wipe out with a dry soft cloth.

TECHNICAL SPECIFICATIONS

Power: 70 W

Supply: 230 V 50 Hz

Noise L_{WA}: 65 dB

KB MAX: 20 min.

NOTE! MPM agd S.A. reserves the right to introduce technical modifications!

PROPER DISPOSAL OF THE PRODUCT (waste electrical and electronic equipment)

Poland This symbol on the product indicates that the product, after its lifetime, should not be disposed with other household wastes. To avoid harmful influence on the environment and human health due to uncontrolled waste disposal, dispose the waste device to the service point of waste household appliances or report collecting it from home. In order to obtain detailed information about the location and how to dispose the waste electrical and electronic equipment in a safe way, contact your retailer or local Department of Environmental Protection. Do not dispose the device with other municipal waste.

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