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
## Product details
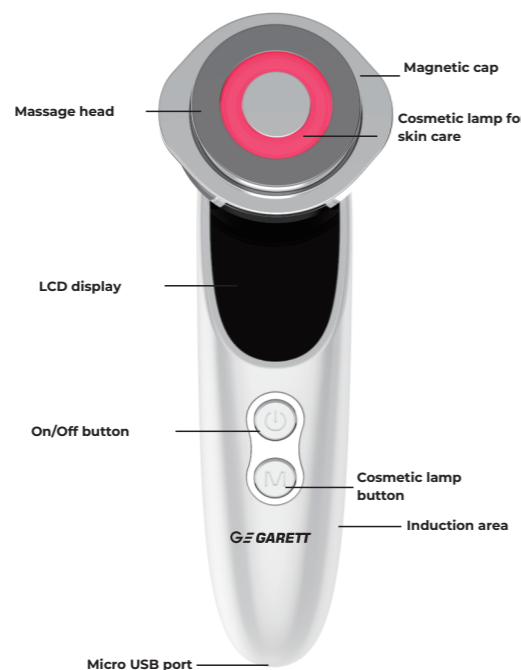
## Function details

**HEATING MODE - (HOT)**  
The heating mode (with vibrations off) cleanses the pores and thoroughly removes make-up from the skin.

**Steps:**

- remove the magnetic cap from the device and place a cotton pad soaked in make-up remover underneath it,
- press the On/Off button to select the HOT mode on the display.
- hold the induction area on the back of the device to activate it and cleanse your face with the device head.

**VIBRATION MODE WITH ELECTRICAL STIMULATION - (EMS)**  
The vibration mode with electrical stimulation is used to tighten the skin, increase its elasticity and improve the metabolism of skin cells.

**Steps:**

- before the treatment, cleanse your face with a toner, essence or skin lotion,
- press the On/Off button to select the EMS mode on the display. Select the appropriate operating mode for the cosmetic lamp with the M button.
- hold the induction area on the back of the device to activate it and begin the treatment.
- move the vibrating head across the face in an outward and upward motion for the desired effect.

**VIBRATION + HEATING MODE - (EMS+HOT)**  
The vibration + heating mode combines both previously described operating modes.

It is perfect for skin care, used alongside skin care products.

**Steps:**


- before the treatment, cleanse your face and apply your favourite skin care product.
- press the On/Off button to select the EMS + HOT mode on the display. Select the appropriate operating mode for the cosmetic lamp with the M button.
- hold the induction area on the back of the device to activate it and begin the treatment.
- move the vibrating head across the face in an outward and upward motion for the desired effect.

**COSMETIC LAMP FOR SKIN CARE (PHOTON THERAPY)**  
The cosmetic lamp also has 3 operating modes:

- R (Red) - whitening and rejuvenating the skin, shrinking enlarged pores and treat broken blood vessels;
- G (Green) - brightening the skin;
- B (Blue) - eliminating acne, reducing sebum production, killing acne-causing bacteria and improving oily skin.

**NOTE:** Remember to put your hand on the induction area on the back of the device to activate the operating modes. In order to achieve the desired results, depending on the selected mode, use a cotton pad, facial cleanser, and moisturising cream or oil.

## LCD display

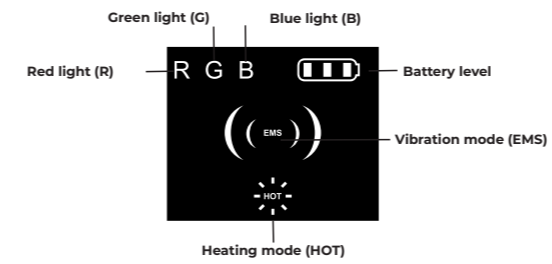**CONTROL PANEL:**  
1. Pressing and holding the On/Off button for about 2 seconds turns on the device in EMS vibration mode.  
2. A short press on the On/Off button activates the HOT + EMS mode.  
3. A short press on the On/Off button activates the HOT mode.  
4. Pressing the On/Off button again sets the mode back to EMS (3 modes: EMS, HOT+EMS and HOT mode).  
5. Press the M button to turn on the skin care lamp. Each successive press of the button changes the colour of the lamp.  
6. The device starts to vibrate when the vibration head comes into contact with the skin. Remember to keep your hand in contact with the induction area on the back of the device.  
7. Press and hold the On/Off button for 2 seconds to switch off the device.

**NOTE:** In each new device, the induction area is protected with film. Until the film is removed from the area, the device will not activate any of the above-mentioned operating modes.

## Product specifications

Product dimensions	167*49*40mm
Power	2.5W
Charging voltage	USB5V
Charging time	90min
Charging capacity	450mAh
Net weight	106g
Set	user manual, micro USB charging cable

## Precautions

1. Do not use the device in case of:
  - previous facelift treatments
  - metal implants under the skin
  - implantation
  - pregnancy
  - menstruation
2. Depending on the selected operating mode, use a cotton pad, facial cleanser, and moisturising face cream or oil when using the device.
3. After each use, thoroughly clean the head of the device. Wash the dirt off the surface with a facial cleanser.
4. Do not use the product on dry skin.
5. The product can only be used on the face and neck.
6. Avoid using the device in case of wounds or skin diseases.
7. People with sensitive skin are recommended to use the device only once a week.

## WARRANTY/RETURNS

The product has a 24-month manufacturer's warranty. The warranty does not cover defects and damages caused by the user's fault or as a result of improper use (incompatible with the intended use). The warranty also does not cover mechanical damage and damage caused by using accessories other than those included in the kit. The warranty also does not cover consumables. When advertising the product, it is necessary to fill in the warranty card attached to the device and send it together with the proof of purchase to the service address given on the warranty card.

Contact details are available on the website: [www.garett.eu](http://www.garett.eu)

**The producer:**  
PP S.C. Paweł Knez, Piotr Wilczyński  
Targowa 18/1413  
25-520 Kielce  
[www.garett.eu](http://www.garett.eu)

A marking on the product or in the texts referring to it indicates that it should not be disposed of with other household waste after its useful life has expired. To avoid harmful effects on the environment and human health from uncontrolled waste disposal, please separate the product from other types of waste and recycle responsibly in order to promote the reuse of material resources as a sustainable practice. For information on where and how to recycle this product in an environmentally sound manner, household users should contact the retailer where they purchased the product or the local authority. This product should not be disposed of with other commercial waste.

Made in P.R.C.